---

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**OVE**

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---

Lucibel·le

PARIS

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THE POWER OF LIGHT

---

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THE OVE BEAUTY MASK BY LUCIBELLE PARIS IS BORN OUT OF A REVOLUTIONARY COMBINATION OF CUTTING-EDGE TECHNOLOGY AND DESIGN, DELIVERING UNRIVALLED BEAUTY.

The result of an unprecedented partnership with designer Olivier Lapidus, its light pioneers a new path to beauty. The technology of the **OVE** beauty mask capitalizes on skin regeneration through LED photobiomodulation.

Red LED light, a source of beauty and wellness, boosts the synthesis of collagen, elastin and hyaluronic acid, giving the skin a replenishing, smoothing and invigorating effect.

To benefit fully from the **OVE** mask's features and to avoid any damage linked to improper usage, carefully read the instructions prior to use.

**A new era :  
beauty through  
light**

17

**P19** SAFETY INSTRUCTIONS

**P20** CONTRAINDICATIONS

**P22** HANDLING THE  
BEAUTY MASK

**P24** OPERATING THE  
BEAUTY MASK

**P25** DIRECTIONS FOR USE  
AND MAINTENANCE

**P26** TECHNICAL  
SPECIFICATIONS

# The power of light

## CONTENTS OF THE BOX

- ✓ An **OVE** beauty mask
- ✓ A pair of protective glasses
- ✓ A user guide

## SAFETY INSTRUCTIONS

To avoid any risk of electrocution, fire, electric shock or injury, please follow the following safety measures:

# 19

**WARNING:** Do not look directly into LEDs while in operation. You must wear the goggles provided, due to the intense optical radiation.

**Check** that the network voltage is in accordance with the range indicated on the power terminal of the **OVE** beauty mask (100-240 VAC).

**Plug** the mask into a wall socket which meets electrical specifications.

**Do not use** the **OVE** beauty mask under running water, while taking a bath or in the shower.

**Never use** the beauty mask with a damaged cord and/or power supply.

**Do not use** the beauty mask if the timer is faulty or if the inner pane of the mask is broken or has been removed.

**Do not place** the mask on a damp surface or close to a bathtub, shower or basin.

**Use** the **OVE** beauty mask in a dry place and avoid any sprayed water.

If the **OVE** beauty mask falls in water: **do not touch it** and immediately disconnect the power.

**Do not use** the beauty mask under a cover or in a similar environment. Excessive heat could be emitted and cause a fire, an electric shock or bodily harm.

**Only carry and move** the **OVE** beauty mask by holding the base. Never lift by the flexible cable.

**Do not shake** or drop the **OVE** beauty mask.

**Do not leave** the **OVE** beauty mask unattended when it is turned on and plugged in.

**Do not remove** the serial number located beneath the mask base.

**Do not open** or dismantle the **OVE** beauty mask. Maintenance must only be carried out by the manufacturer in order to avoid any danger.

**In the case** of use of the mask by a third party, the instructions for use must be provided.

**WARNING :** Do not stare at the light sources.

It is necessary to wear the glasses provided because of the light intensity. These goggles must be worn during the entire duration of the photobiomodulation sessions.

Cleaning and user maintenance shall not be made by children without supervision.

The **OVE** beauty mask is designed for use on the face only. In the event of improper handling, in particular for use of the mask on other part of the body, the legal guarantee of conformity, which the consumer enjoys, shall not apply.

## CONTRAINDICATIONS

The **OVE** beauty mask is suitable for all skin types.

**However, its use is not recommended for:**

- Pregnant women,
- People suffering from telangiectasia and photo-dermatosis,
- People with birth marks,
- People who suffer from or who may suffer from epileptic fits,

Warning! Not to be used by persons under the age of 18 years.

**This appliance** must not be used by individuals with impaired physical, sensory or mental capacities, or those lacking experience and knowledge of the appliance's use and its potential risks.

Children should be supervised to ensure that they do not play with the appliance.

Children shall not play with the appliance.

In addition, use of the **OVE** beauty mask is not recommended for people undergoing treatment with photosensitizing medications or people being treated for a disease which results in photo-sensitivity.

Pre-existing skin conditions cannot be attributed to exposure to the mask's red light.


The use of the OVE face mask is contraindicated for persons suffering from sunburn.

If undesirable side effects are observed, such as itching within 48 hours of initial use, please ask your doctor prior to reusing the appliance.

Use of red light is not contraindicated with anti-cancer treatments.

If in doubt, user should ask a doctor.



# 21

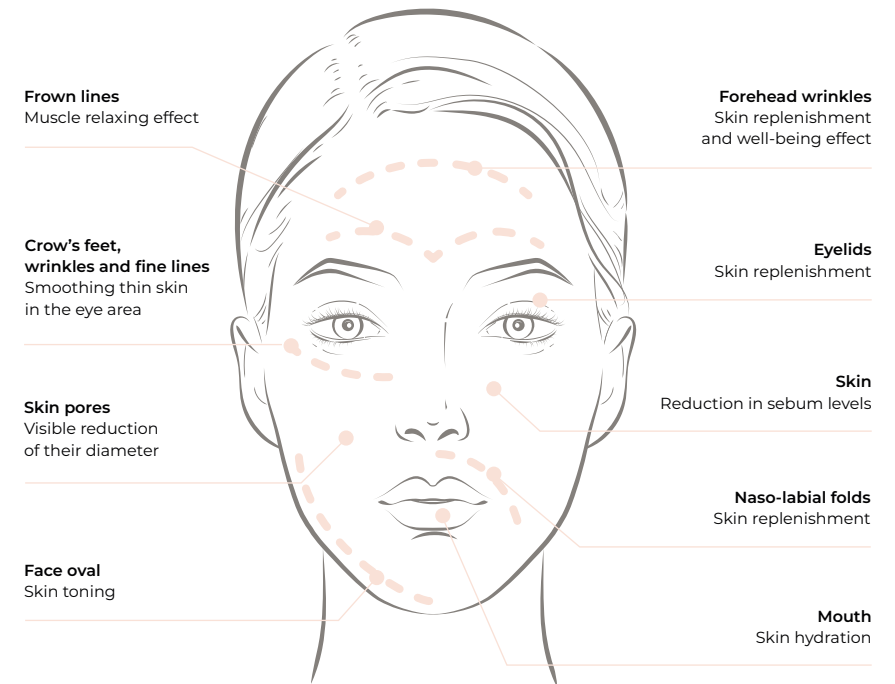
**Thanks to the LED light intensity, specific to every area of the face,** the beauty mask offers unparalleled performance for visible results.

# HANDLING THE BEAUTY MASK


Advice for simple and efficient use of the OVE beauty mask.

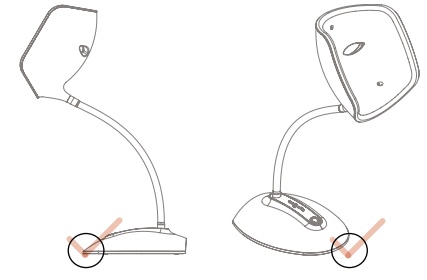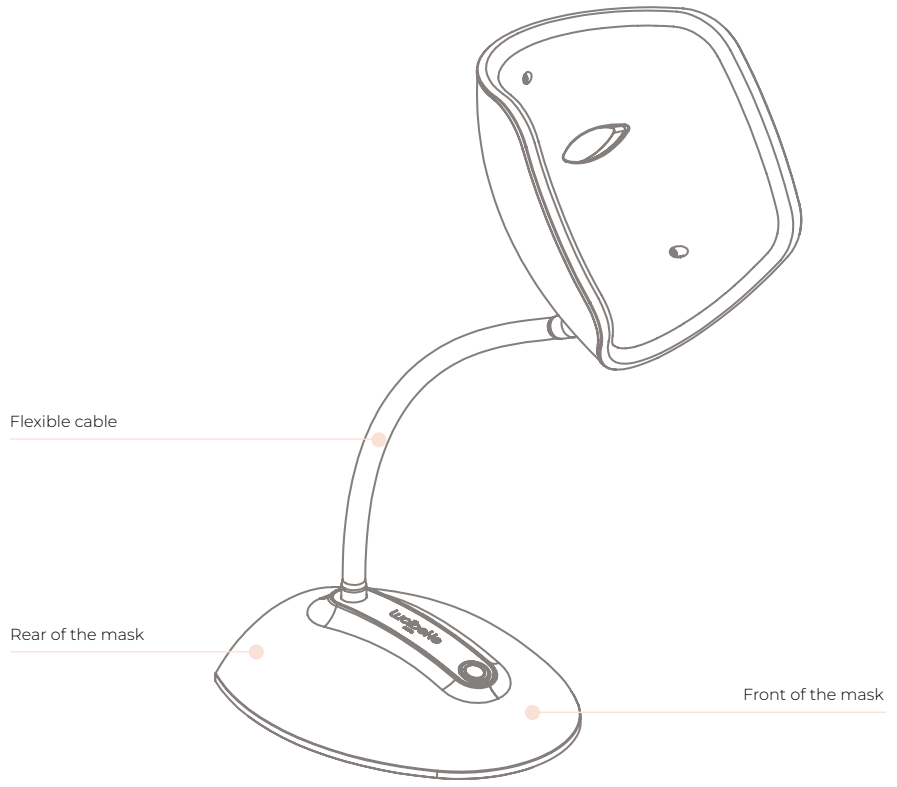
Always turn the mask the right way around in relation to the base. The light emitting area of the beauty mask must always be facing the front of the base. To keep the mask facing the right way with regard to the base, rotate the flexible cable.

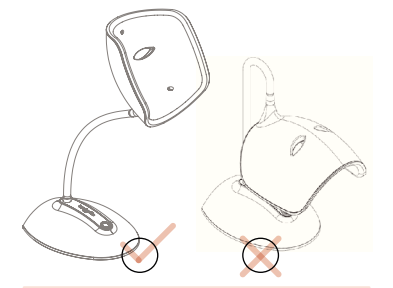Ensure not to over-bend the cable, as you may risk damaging or knocking over the mask.

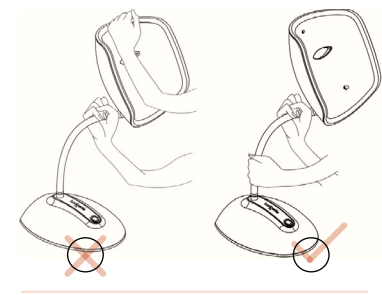To adjust the height of the mask, always handle the flexible cable using both hands.

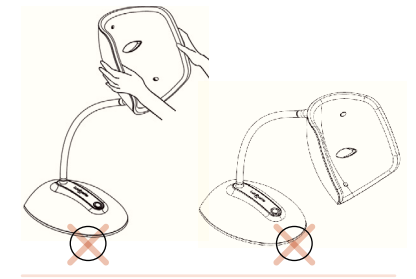The beauty mask can only be used in one position. Do not rotate: adjust the mask using the flexible cable.

## OPERATING THE BEAUTY MASK

**Place** the mask on a flat, stable surface and plug it into a power outlet (for optimal and safe use, we do not recommend using a socket strip).

**Adjust** the protective goggles on the face so that the eyes are fully protected during photobiomodulation sessions.

**Place** the face inside the mask. The nose must touch the wall. If needed, adjust the position of the mask by moving the flexible cable.

**Press** the button once (short press) to turn the beauty mask on alone.



The session is programmed to last for **12 minutes**, alternating flashes and steady light.

To use a Lucibel.le Paris accessory alone: plug it into the jack located behind the mask and press the button once (long press).

To use an accessory and the mask at the same time, press the button twice.

After the session, the mask turns off automatically. At the end of the session, the mask and the accessory are automatically turned off. To stop the session, press the button briefly.

**After use**, disconnect the mask from the power supply.
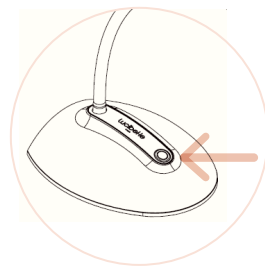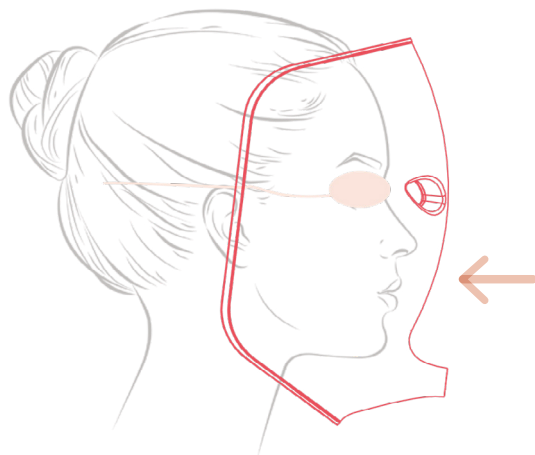
## DIRECTIONS FOR USE

**Cleanse the skin** prior to use to optimize light absorption.

Set 72-hours **intervals between sessions** - the time needed for the skin to absorb the light energy.

Do not use the **OVE** face mask for longer than the recommended time.

Optimal effects are observed after 24 sessions on average. To benefit from the “**delayed effect**”, take a 1-2 month break before resuming sessions. During this pause time, the skin will continue to be embellished.

**Only use** accessories supplied by Lucibel.le Paris.

Avoid applying any cosmetic products before the session other than those supplied by Lucibel.le Paris.

Finally, for an accurate self-assessment of the beauty mask's results, take a photo of the treated areas prior to starting the protocol, a second photo after 16 sessions and a third photo after 24 sessions.

If the face mask is used by several people in succession, keep it off for at least 5 minutes between applications.

The mask is designed for home use with up to five sessions per day; more intensive use would damage the beauty mask.

# 25

## Maintaining the beauty mask

Clean the beauty mask after each use using wipes soaked in an antibacterial and antiviral solution.

Refrain from using detergents, solvents, chemical products and abrasive sponges.

## Environmental safeguards

The beauty mask may not be treated as household waste: it must be brought to a designated collection point for the recycling of waste electric and electronic equipment, to help preserve natural resources.

## TECHNICAL SPECIFICATIONS

Technical specifications of the beauty mask - a real marvel of technology - exalting every kind of beauty.

<b>Operating temperature</b>	0 °C to +35 °C +32 °F to +95 °F
<b>Maximum service life</b>	50,000 hours
<b>Materials</b>	Casing: PC, ABS Flexible cable: PVC, Brass Base: Steel Screws and bolts: Zinc-plated steel Dissipator: Aluminum
<b>Weight</b>	3.64 kg
<b>Dimensions</b>	W 213 mm x L 396 mm x H 611 mm
<b>Number of LEDs</b>	40
<b>Electrical input characteristics</b>	100/240 V ~ 50/60 Hz
<b>Wavelength</b>	635 nm
<b>Power</b>	40 W
<b>Jack socket</b>	Reserved for Lucibel.le Paris accessories

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